

Plate Count Agar (Standard Methods Agar)

Intended Use

Plate Count Agar (Standard Methods Agar) is used for obtaining microbial plate counts from milk and dairy products, foods, water and other materials of sanitary importance.

Summary

Plate Count Agar is formulated as described by Buchbinder *et al.*, which is recommended by APHA and FDA.

Principle

Tryptone provides nitrogenous and carbonaceous compounds, long chain amino acids, and other essential nutrients. Yeast extract supplies Vitamin B complex. APHA recommends the use of pour plate technique. The samples are diluted and appropriate dilutions are added in petriplates. Sterile molten agar is added to these plates and plates are rotated gently to ensure uniform mixing of the sample with agar. The poured plate count method is preferred to the surface inoculation method, since it gives higher results. Plate Count Agar is also suitable for enumerating bacterial count of sterile rooms.

Formula*

Ingredients	g/L
Tryptone	5.0
Yeast Extract	2.5
Dextrose	1.0
Agar	15.0
Final pH (at 25°C)	7.0 ± 0.2

*Adjusted to suit performance parameters.

Storage and Stability

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Avoid freezing and overheating. Use before expiry date on the label. Once opened keep powdered medium closed to avoid hydration.

Type of Specimen

Food, Milk and Dairy samples; Water samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines.

For food and dairy samples, follow appropriate techniques for handling specimens as per established guidelines.

For water samples, follow appropriate techniques for handling specimens as per established guidelines and local standards.

Specimens should be obtained before antimicrobial agents have been administered.

After use, contaminated materials must be sterilized by autoclaving before discarding.

Directions

1. Suspend 23.50 g of the powder in 1000 mL purified / distilled water.
2. Mix thoroughly.
3. Boil with frequent agitation to dissolve the powder completely.
4. Sterilize by autoclaving at 121°C (15 psi) for 15 minutes as per validated cycle.

Quality Control

Dehydrated Appearance: Light yellow coloured, homogenous, free flowing powder.

Prepared Appearance: Light yellow coloured, clear to slightly opalescent gel forms in petridishes.

Growth Promotion Test: Growth promotion is carried out in accordance with the harmonized method of USP/EP/JP/IP and growth is observed after an incubation at 30°C-35°C for 18 - 48 hours.

Growth Promoting Properties: The test results observed are within the specified temperature and shortest period of time specified in the test, inoculating ≤100 cfu of appropriate microorganism at 30°C-35°C for 18 hours.

Organism (ATCC)	Growth
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> (6538)	Good
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> (25923)	Good
<i>Escherichia coli</i> (8739)	Good
<i>Escherichia coli</i> (25922)	Good
<i>Pseudomonas aeruginosa</i> (9027)	Good
<i>Bacillus spizizenii</i> (6633)	Good
<i>Enterococcus faecalis</i> (29212)	Good
<i>Streptococcus pyogenes</i> Strain Bruno (19615)	Good

Note: For good growth - Growth obtained on test media should not differ by a factor greater than 2 from calculated value for a standardized inoculum.

Performance and Evaluation

Performance of the product is dependent on following parameters as per product label claim:

1. Directions
2. Storage
3. Expiry

Warranty



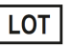


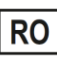



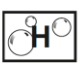
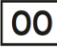
This product is designed to perform as described on the label and package insert. The manufacturer disclaims any implied warranty of use and sale for any other purpose.

Reference

1. Buchbinder L., Baris Y., Aldd E., Reynolds E., Dilon E., Pessin V., Pincas L. and Strauss A., 1951, Publ. Hlth. Rep., 66:327.
2. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
3. FDA Bacteriological Analytical Manual, 2005, 18th Ed., AOAC, Washington, DC
4. Data on file: Microxpress®, A Division of Tulip Diagnostics (P) Ltd.

Product Presentation:

Cat No.	Product description	Pack Size
201160210100	Dehydrated Culture Media	100 g
201160210500	Dehydrated Culture Media	500 g
203160470100	Bottle Media	100 mL
203160470250	Bottle Media	6 x 250 mL
203160470500	Bottle Media	500 mL
205160510100	Ready Prepared Plate (90mm)	100 plates

 Temperature Limit	 Manufacturer	 Batch Code	 Date of Manufacture	 This way up	 Received on
 Catalogue Number	 Consult Instructions for use	 Use-by Date	 Hygroscopic keep container tightly closed	 Opened on	

Revision: 0725/VER-03

Disclaimer

Information provided is based on our inhouse technical data on file, it is recommended that user should validate at his end for suitable use of the product.